

# SOLID PLUS HEAT

[WWW.ROTAID.COM](http://WWW.ROTAID.COM)



**TELEFUNKEN**

 **ROTAID™**  
AED CABINETS

# ROTAID SOLID PLUS HEAT

## PRODUCT SHEET

The Rotaid Solid Plus Heat is an innovative AED wall cabinet that ensures a recognizable and robust storage of your AED device. The characteristic round shape and green colour, make it easy to be found and recognized and comes with a high quality appearance. The large AED logo, positioned in the transparent centre of the cover, leaves no room for interpretation. This is a life saving device.

The cabinet can be opened by an intuitive twist that removes the whole cover. The responder therefore has full access to the AED that hangs via a bracket for a quick and easy grab. The cover is attached to the back plate via a belt and needs no further attention from the responder. A stainless steel hook, attached to the bottom of the cabinet, can hold a seal, tag or other appliance.

The cover, equipped with recesses for an optimal grip and directive arrows, is made from 100% Polycarbonate. With a thickness of 5 millimetres it offers a solid protection and ensures a high grade of water- and dust resistance. Furthermore, the cover is fitted with UV filters, to protect from direct and indirect sunrays. Appearance decay, rust and accumulation of residue are excluded.



A tailor made heating device maintains an inside temperature of at least 5 Degrees Centigrade even when temperatures drop to -25 Degrees Centigrade. The automatic and gradual operation makes it particularly effective and efficient.

The materials are tested in the outdoors and provide a long lasting protection of the AED.



**VERY LOW  
POWER  
CONSUMPTION**


**FULLY AUTOMATIC  
OPERATION**

**OPERATIONAL AT  
TEMPERATURES  
AS LOW AS  
-25°C**

**LED LIGHTING**

# ROTAID SOLID PLUS HEAT

## TECHNICAL DETAILS

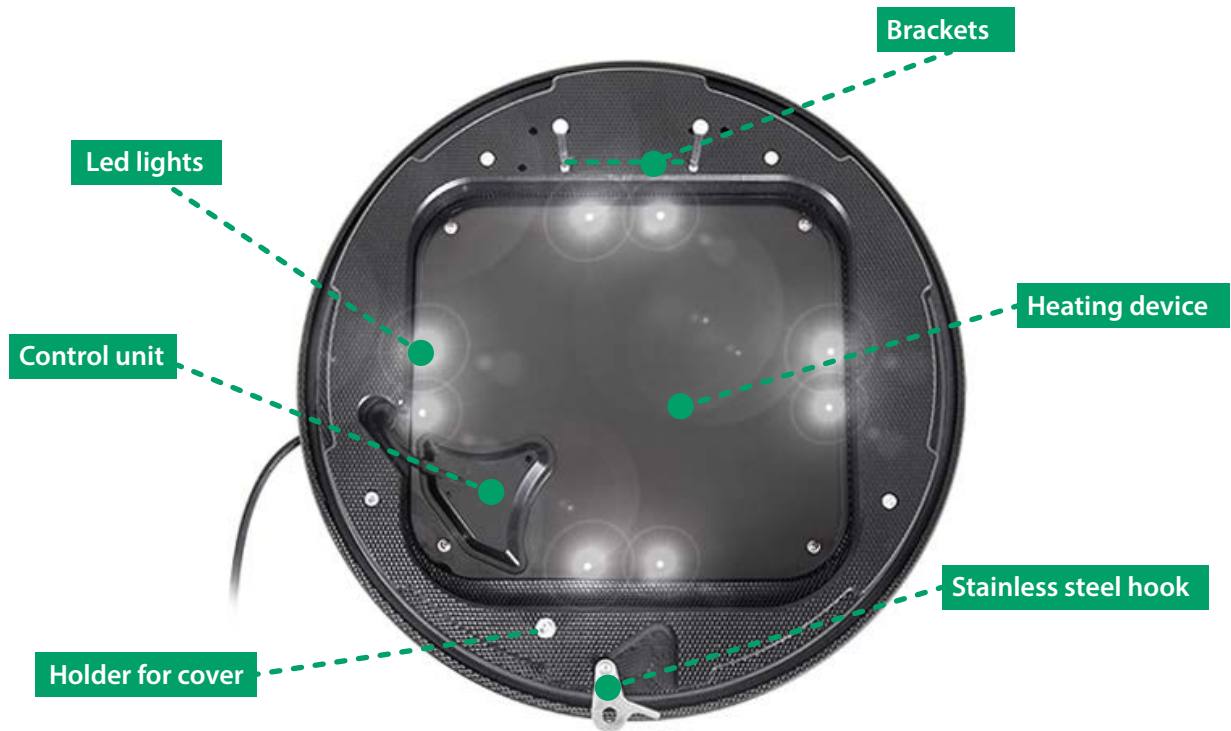

The back plate is made from 100% HDPE plastics and offers space and protection to practically all types of AED devices on the market. An optional audible alarm (90dB) is available.



### Dimensions

Length **48 cm** Width **48 cm** Depth **18 cm** Weight **6,7 kg**

Length **18,9"** Width **18,9"** Depth **7,1"** Weight **14,7 lbs**

Packed individually by box, fastening materials included.

[www.rotaid.com](http://www.rotaid.com)